promed pedisenso duo

Promed GmbH
Lindenweg 11
D-82490 Farchant
Tel: +49 (0)8821/9621-0
Fax: +49 (0)8821/9621-21
info@promed.de
www.promed.de

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SE

FI

- A** Mains adapter
- B** Interlock release
- C** Power plug
- D** Device
- E** Collet chuck adapter support
- F** Speed adjustment / On-Off switch
- G** Collet chuck adapter
 - 1. Safety collet chuck adapter for diabetic (yellow)
 - 2. Collet chuck adapter (white)

Please keep this user manual in a safe place!

Please read the user manual thoroughly before using the appliance for the first time.

Thank you for purchasing the **Promed pedisenso duo** manicure and pedicure appliance.

You have opted for a German quality product of a leading manufacturer of skin and nail care appliances. Expert advice is available at all places where you can buy **Promed** products or feel free to contact us with your queries. We can recommend a consultant. The equipment was developed specially for the use as manicure and pedicure appliances in the domestic private use.

DESCRIPTION OF THE UNIT!

The Promed pedisenso duo is specially suitable for a diabetic due to its power transmission with a planetary gear and its yellow safety collet chuck adapter. When you purchase the **Promed pedisenso duo** you get a top-quality and elegant manicure/pedicure set, with luxurious ergonomics and professional-quality filing accessories. It is easy to use, making it ideal for hand, foot and nail care, especially for ingrowing, hardened and thickened nails and for treating corns, calluses and hard skin. It is also suitable for cosmetic care for your hands, fingernails and artificial nails. The promed emotion set offers gentle and reliable hand and nail care and is extremely practical to use. The device is equipped with 2 collet chuck adapters which are provided for numerous cutting attachments and accessories in dental measurements (Ø 2.332 - 2.350 mm).

We, as the manufacturer, cannot be made liable in any way for injury or damages to people or objects that arise from failure to comply with this warning.

We wish you a lot of enjoyment with your new **Promed pedisenso duo** unit. We would like to familiarize you with the unit in the following sections. Please read the user manual thoroughly before using the unit for the first time.

Scope of delivery:

Device including cable, transformer plug, case, cleaning brushes, 2 collet chuck adapter and a 41-piece grinding set.

Additional attachments are also available from specialist shops. Inform yourself regarding the details.

DE

EN

FR

IT

ES

NL

RU

PL

SE

FI

INFORMATION ON THE UNIT / ACCESSORIES

DE

EN

FR

IT

ES

NL

RU

PL

SE

FI

Before operating the unit for the first time, please check that all contents have been delivered and that the line voltage is correct. The supplied power adapter **(A)** is an electronic power supply unit intended for an operating voltage of 100 V~ to 240 V~. Should you require a different plug type for the operation of the device, you can purchase it from **Promed**. To change the plug type, press the release **(B)** and pull the power plug **(C)** from the power supply unit **(A)**. Take the required power plug **(C)** and push it back into the power supply unit **(A)**.

- The **Promed pedisenso duo** unit uses a high performance motor. The electronic speed control is used to smoothly set the speed of the grinding head.
- The collect chuck adapter accept all conventional abrasives with shaft in standardised dental measurements of 2.332 mm to 2.350 mm. To ensure a long life span of the collet chuck adapter, please use only quality abrasives with rounded shaft closure and a maximum shaft length of 36 mm!

- Always make sure that the grinding head does not come in contact with other objects when you lay the handpiece down. There is a risk of the shaft becoming compressed and asymmetrical. Asymmetry in the grinding head causes damage to the handpiece!
- Grinding head bits are available at your local specialist retailer or **Promed**. (www.promed.de).
- To ensure hygienic operation, the grinding heads must be disinfected after use.
- The handpiece may only be disinfected using disinfectant wipes.

INFORMATION ON THE UNIT / ACCESSORIES

- Information about suitable disinfectants and products can be obtained from **Promed**.

For disinfection of instruments and grinding heads:

Promed Pure ID (Product No. 330810)

For disinfection of surfaces and devices:

Promed Pure FD (Product No. 330815)

For disinfection and cleaning of instruments and grinding heads: **Promed UC-50** (Product No. 330210)

OPERATING THE UNIT

Before you switch on the device, you must regulate the speed (rotary knob **(F)**) to zero.

The power plug is intended as a separator. Make sure that the power plug is plugged into an easily reachable wall socket.

- Plug the power adapter **(A)** into a power outlet.
- Select a collet chuck adapter **(G)** and screw it into the collet chuck adapter support **(E)**.
- Select a suitable grinding head. The yellow safety collet chuck adapter **(G)** is a safety collet chuck adapter with which you may plug the abrasive at the front side of the adapter. Once the abrasive body begins to rotate, the mechanism locks automatically. With the white collet chuck adapter **(G)** you must open the knurled screw by turning it, then plug the abrasive and lock the chuck by tightening the knurled screw again. Never use force to remove or insert a grinding head!

DE

EN

FR

IT

ES

NL

RU

PL

SE

FI

OPERATING THE UNIT

Abrasives with 2.332 mm to 2.350 mm shank diameter fits precisely with the collet chuck adapter with ceramic bond, without you having used excessive force. Please switch off the appliance for removing an abrasive. Now remove the abrasives by simply pulling out straight. Now remove (at the white adapter after loosening the knurled screw) the abrasive by simply pulling it out.

D. Speed adjustment and On-Off switch

On the rear part of the device, you will find the speed control **(F)**, with which the grinding force can be regulated smoothly. At the same time the speed control **(F)** serves as on-off switch of the device.

E. Before initial application, run the hand piece necessarily for approx. 30 minutes with approx. 2.000 rpm (approx. ½ revolutions of the controller). This is necessary to bed in the carbon brushes of the motor.

Note:

A different running noise between clockwise and counter-clockwise rotation is normal and caused by bedding in of the carbon brushes or by running in of bearings.

In cases of malfunctioning:

- Check that the power supply is correct.
- Check if the collect chuck adapter is screwed properly in the collect chuck adapter support.
- If the problems have still not been solved after following these steps, please consult your service center immediately.

DE

EN

FR

IT

ES

NL

RU

PL

SE

FI

SPECIAL QUALITY FEATURES

Planetary gear unit:

Through the planetary gear, the device can be transferred at a low engine speed with much more effective force and torque on the abrasive. This makes the device specially suitable for a diabetic, as a result the risk of injury is reduced to minimum.

Safety collet chuck adapter:

The yellow safety collet chuck adapter serves a unique safety system that stops the revolutions of the abrasive with much pressure

Safety Stop:

If the device develops excessive speed and pressure during grinding, the electronics will respond by switching off the device. In this case, you must set the speed control **(F)** to zero to fix the device back into operation.

USAGE

Skin and nail care is not merely a matter that is indicated when suffering correspondingly, it is also a cosmetic matter.

Start with cutting and forming toenails. They absolutely must be filed straight in order to hinder nails growing into the skin. Micro-fine filing using the sapphire disc protects against ingrown nails. The micro-fine files with the sapphire disc prevent an ingrowth. Troublesome calluses always develop in the areas where the foot is over-stressed, mostly the heel, the front balls and the side of the big toe. Treat the callus using the sapphire cone or even better with the drainage grinder. Repeat this procedure regularly. Under no circumstances should you take radical measures when removing the callus.

Ingrown toenails can – as needed – be loosened gently using the small cutter. Thick lignified nails can be ground off slowly from the top and filed thinner from the sides.

DE

EN

FR

IT

ES

NL

RU

PL

SE

FI

USAGE

DE

EN

FR

IT

ES

NL

RU

PL

SE

FI

Crumbly nail areas can also be cleaned using the fine sapphire cutter.

Corns are the painful result of a slow callus process. They are calluses with an externally hardened part which is soft inside. Remove the upper callus with the sapphire cone and then work further using the small sapphire cutter. This already removes the pressure. Deeply sunken corns should not be removed by oneself; they should be removed by a doctor.

If skin and nail care are carried out systematically using the **Promed** devices and accessories, regular care becomes a positive experience.

The special abrasive wheels are described more closely on the following pages.

The abrasive wheels are available in various forms and degrees of hardness and can be reordered at any time. Ask your specialist dealer.

Note:

Do not wash your hands or, when applicable, your feet before treatment. A soothing and calming moisturizer should be first used after treatment.

APPLICATION OF THE ABRASIVE WHEELS

Stone sander, fine, pointed

For the removal of difficult corners and edges at those sites that thicker grinding heads cannot access. (Be careful when using the needle burr!) To clean the underside of fake nails.

Stone sander, cylindrical

For the grinding of tough, thickened nails. Thickened nails can be shaped easily and their appearance improved. Pressure is relieved. It is also suitable for grinding the edges of hard corns to relieve pressure. (Only grind away the areas of pressure!)

Sapphire needle burr, pointed

For the removal of difficult corners and edges that cannot be accessed by thicker grinding heads. (Be careful when using the needle burr!)

Sapphire file, cylindrical

For the grinding of tough, thickened nails. Thickened nails can be shaped easily and their appearance improved. Pressure is relieved. It is also suitable for grinding the edges of hard corns to relieve pressure. (Only grind away the areas of pressure!)

DE

EN

FR

IT

ES

NL

RU

PL

SE

FI

APPLICATION OF THE ABRASIVE WHEELS

DE

EN

FR

IT

ES

NL

RU

PL

SE

FI

Milling cutter sapphire bit, spike cone rounded

To grind nails infected with fungus, it is very good to remove nail fungus traces. To remove cuticle.

Milling cutter sapphire bit, small spike cone

To grind nails infected with fungus, it is very good to remove nail fungus traces. To remove cuticle.

Callus drainage grinder, sapphire

Particularly suitable for tough, thick calluses on the heel and the ball of the foot, the sites where thick calluses are often present.

Grinding cone, sapphire

For callused and rough skin on the feet and hands, but also for abrading very thick, hard nails. It reshapes crooked nails and relieves painful pressure.

APPLICATION OF THE ABRASIVE WHEELS

Grinding mandrel, plastic + abrasive band, medium

As preparation for the application of fake nails using gel, acrylic, etc.


Grinding disc, sapphire

Shortens and shapes nails right into the corners. Result: perfectly cared for nails that no longer tear, break off or split.

Polishing cone, felt

The polishing cone for the final finish in every type of nail treatment. Toenails and fingernails are buffed to a natural shine. Smooths the edges and surfaces of the nails for the polished finish to the treatment.

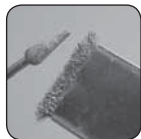Cleaning brush

To clean the grinding heads apart from the abrasive band.

DE

EN

FR

IT

ES

NL

RU

PL

SE

FI

MAINTENANCE AND CARE

Cleaning the unit

- Regular maintenance and cleaning ensure optimal results and extends the life of your unit.
- Store the **Promed pedisensoduo** in a safe, cool and dry location when it is not in use.
- To ensure that the safety, reliability and performance of the unit are not adversely affected, enhancements, readjustments, alterations or repairs are not permitted aside from those performed by **Promed** authorized persons.
- Before cleaning the unit, make sure that it is turned off and disconnected from the electricity.
- Use a slightly moist, soft cloth to clean the unit.
- Do not use any aggressive substances such as scouring agents or other detergents as they may corrode the surface.
- Wipe clean with a soft, moist cloth and do not use scouring agents.
- Do not immerse the unit in fluid and keep the unit away from solvents and strong detergents.
- Do not repair the unit. The unit does not have parts that consumers can service.

Cleaning the grinding heads

The abrasives can be cleaned from dust and skin residuals by using a cleaning brush (e.g. **Promed cleaning brush** Art.No. 198910). It is always advisable to disinfect the abrasives with a commercially available instrument disinfection (e.g. **Promed PURE-ID**), in order to avoid transmission of infections. Please make sure to keep the pins of the abrasives free of grease and clean these from time to time.

DE

EN

FR

IT

ES

NL

RU

PL

SE

FI

SAFETY INSTRUCTIONS

Please pay attention to the following before first time operation:

Never place the unit in direct contact with fire, gas, oxygen or hot objects such as hot plates.

Please observe the following during operation:

A. Grinding heads to be used:

- Standardized dental size: 2.332 - 2.350 mm / 0.091 - 0.093 inches.
- Only use centered grinding heads with a rounded shaft end.
- The use of rusted or bent grinding heads is not permitted.

B. Use the speed controller carefully.

C. Switch the unit off before you change the direction of rotation.

D. Appliance power supply: 100 V~ till 240 V~.

E. Switch the unit off before you change the grinding head.

F. Take every precaution possible to ensure that the unit does not fall down or get damaged in any other way.

G. If problems do arise with the handpiece, take it immediately to get it repaired.

H. Do not lubricate or wash the handpiece.

I. Only for home use.

J. Scrupulous supervision is required when the device is being used in the vicinity of children.

DE

EN

FR

IT

ES

NL

RU

PL

SE

FI

SAFETY INSTRUCTIONS

DANGER!

- Do not place the unit in direct contact with water or other fluids.
-  Never place or use the unit in a wet or moist environment.
- Do not touch the unit with wet hands.
- Operate the appliance only with an appropriate power adapter.
- It is not permitted to use the unit outdoors.

WARNING!

- Make sure that children do not play with the device.
- Never leave the unit unattended when children or inexperienced persons are in its vicinity.
- Only use the unit for the applications described in this manual.
- Only use attachments and accessories from the manufacturer.
- Never use the unit after the electrical lead has been damaged,

when the unit has fallen down or when it has come in contact with water.

- This unit is not designed to be used by people (including children) who have limited physical, sensory, or intellectual abilities or who do not have sufficient experience and/or sufficient knowledge. This does not apply if these people are supervised by a person responsible for their safety or receive instructions from them on how to use this unit.
- Never let children play with the packaging material; there is a risk of suffocation.
- An electrical appliance is not a toy. Children do not recognize the risks that could develop when using electrical appliances. Therefore, operate and store the unit out of the reach of children.
- Do not place or use the unit on hot surfaces or near naked gas flames to ensure that the casing does not melt.

DE

EN

FR

IT

ES

NL

RU

PL

SE

FI

SAFETY INSTRUCTIONS

- Never let the unit come into contact with water. Use a moist cloth to clean the unit. Do not use solvents.
- Never carry or pull the unit by the power cord. Use the plug to pull the power cord out of the wall socket. Do not pull on the power cord itself.
- Do not pull the power cord over sharp edges, wedge it in between objects, hang it down or position or hang it over hot objects (e.g., hot plates) or naked flames. Protect the power cord from heat or oil.
- Only use extension cords when you have satisfied yourself beforehand that they are not defect or damaged.
- If the power cable of the appliance gets damaged, have it replaced by an expert immediately.
- Do not store the unit near a wash basin or bath tub as there is the danger that the unit may fall or be pulled into the wash basin or the bath tub.

- To avoid all risks of injury, blockage or damage to the unit, make sure that the unit does not come into contact with clothing, cords, shoelaces, bristles, wires, threads, hair etc. when in operation.

Note:

The AC adapter may become slightly warm. This does not indicate a malfunction. The device can be used easily.

DE

EN

FR

IT

ES

NL

RU

PL

SE

FI

DISPOSAL

DE

Power tools, accessories, and packaging should be recycled in an environmentally suitable manner.

EN

Only for EU countries:

Do not dispose of power tools with the domestic waste!

In accordance with the European Directive 2002/96/EC on waste electrical and electronic equipment and its

FR

IT

ES

implementation into national law, nonfunctioning power tools must be collected separately and recycled in an

environmentally suitable manner.

NL

This symbol indicates that disposal of this product in the domestic waste is not permitted within the EU. Waste devices contain valuable recyclable material that should be recycled. Recycling also prevents uncontrolled waste disposal from damaging the environment and human health. Therefore, please dispose of all waste devices in the appropriate collection systems or send the device to the place where you bought the device for disposal. They will then recycle the device.

RU

PL

SE

FI

WARRANTY

This unit has been manufactured with the utmost of care and has been examined in detail before leaving the factory. We therefore provide a warranty of 24 months after the date of purchase when presented with the warranty card issued with this appliance, in accordance with the following conditions.

- Within the warranty period, we will replace all defect parts of the appliance free of charge, including labor costs for repairs under warranty, in cases of verifiable faults in material or manufacture that have occurred with correct usage and have been identified during the warranty period.
- This warranty is void when the unit has been improperly used (e.g., connection to unsuitable sources of electricity, breakage), the unit has been tampered with (e.g., opening the unit's casing), and when spare parts have been used that have not been approved of by **Promed**. Parts subject to wear are exempt from the warranty. Appliance motors and mobile parts are not included in the warranty.
- The warranty period begins on the date of purchase. Warranty

claims do not affect the warranty duration. Warranty claims must be lodged within the warranty period. Claims arising after cessation of the warranty period cannot be considered.

- The warranty only comes into effect within the scope of these warranty conditions when the date of purchase has been confirmed by the stamp/ signature of the dealer on the warranty card.
- In cases of warranty or repair claims, please post the complete appliance with the completed warranty card to the customer service point responsible for you. The warranty card can be found at the end of this user manual.
- We reserve the right to make technical and optical changes as well as changes to the equipment!
- The following information is vital for problem-free processing:
 1. Original sales slip/receipt or dealer's stamp with the date of purchase
 2. Defect detected
 3. Unit name / type

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Symbol for protection class II appliances

The symbol corresponds to symbol 5172 of the EN 60417-1:2000 Standard. The appliance is double insulated as a protective measure. This prevents conductive parts of the appliance from becoming impaired, something that could lead to malfunctioning of the basic insulation voltage.

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Symbol for WARNING! Indicates sections of the user manual that are relevant for safety!

The symbol corresponds to symbol 0434 of the DIN ISO 7000:2008 Standard. Instructions for an appliance, system or accessory that contain important information for the user that are relevant for safety.

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Pay attention to the symbol for the user manual!

The symbol corresponds to the symbol 1641 in the DIN ISO 7000:2008 Standard. It is imperative that the manual be read before using the unit.

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Only transport in the original
packaging.

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EN

Transport- und Lagertemperatur:
+5° C – +50° C

Transport and storage temperature:
+5° C – +50° C
(+41° F to +122° F)

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Relative humidity:
30% to 75% non-condensing

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An abrupt change in temperature is
to be avoided as condensation may
arise in the unit. Only operate the
unit after the temperature has been
equalized.

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EN | The power supply unit complies with the Ecodesign Directive (ErP) No. 2009/125/EC requirements for external power supply units and the requirements of phase 1 and 2 of the EC Directive 2009/125/EC.

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